

9272030

Wireless Gen2 Pull Cord Surface Mount

Description


The Wireless Gen2 Pull Cord Surface Mount is a reliable pull cord solution that can transmit alerts wirelessly to a miPOSITION enabled miALERT Nurse Call System using the proprietary miNET wireless transmission standard. Alerts are initiated using the pull cord and can be cancelled with a physical button the unit face.

The surface mount design allows the unit to be installed directly to the wall with no cutting. It can also be installed on a 1-gang box. Alerts can be transmitted up to 100ft under nominal conditions. Low battery reporting is available through the miALERT Event Management Software.

These pull cords can also be used with Gen2 wireless enabled miALERT voice to room devices.

Features

- Physical cancel from the unit.
- Surface mount design requires no wall cutting.
- Wireless transmission requires no physical cabling.
- Battery powered, no direct power required.
- Two year battery life under nominal conditions.

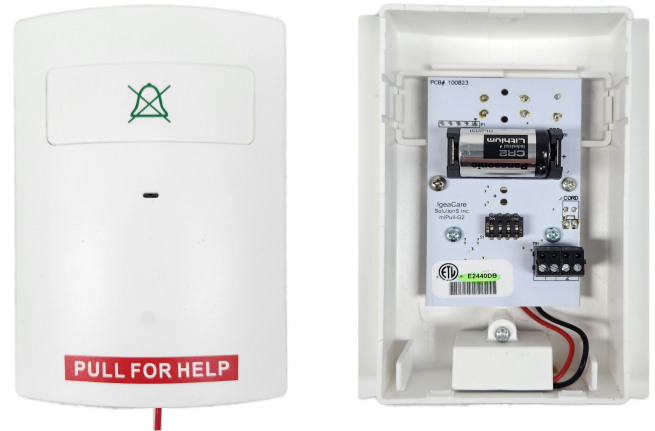

Specifications

Unit with mounting plate: 130mm x 92mm x 33mm

Power: Single 3V CR2 battery cell. Expected battery life under normal conditions is approximately 2 years. Units should be tested at least once a month to ensure battery condition is nominal. For low battery reporting a miALERT Event Management Software server must be present.

Product is designed to conform to the UL 2560 standard for wireless nurse call devices.

For replacement batteries, contact NSI Distribution (www.nsidistribution.com).

Package Contents

- Wireless Gen 2 Surface Mount Pull Cord
- Mounting Plate
- Two (2) 6/32x1 Gang Box Screws
- One (1) 3/4 Retaining Screw
- One (1) CR2 3V Lithium Battery